



## Diaphragm Pumps for Air, Gases and Vapours



### Series LABOPORT® N 811 K\_.18 Pumps

## LABOPORT® Mini Diaphragm Vacuum Pumps and Compressors

#### Technical features:

- 100% oil-free transfer
- Pure transferring, evacuation and compression of gases
- Version for slightly aggressive or corrosive gases and vapours
- Maintenance-free
- Environmentally friendly
- High level of gas tightness.

Series N 811 diaphragm pumps are single-head, dry-running devices used in a wide range of laboratory applications. They transfer, compress and pump down without contamination. The heart of these very compact pumps is a KNF structured diaphragm. This patented diaphragm was stress-optimized using the Finite Elements method. As a result, we were able to make the pumps smaller while increasing the service life of the diaphragm.

The pumps are available in various versions differing in the materials which contact the media.

#### Technical data:

	N 811 KN.18	N 811 KT.18
Delivery (l/min) <sup>1</sup>	11.5	11.5
Ultimate vacuum (mbar abs.)	240	290
Operating pressure (bar g)	2	2
Connectors for tube (mm)	ID 6	ID 6
Permissible gas and ambient temperature	+5...+40 °C	+5...+40 °C
Mains	230V/50Hz	230V/50Hz
Motor protection	IP 20	IP 20
Power P <sub>1</sub>	65 W	65 W
Operating current	0.8 A	0.8 A
Weight	2.5 kg	2.5 kg
Dimensions LxHxW (mm)	187/157/90	187/157/90
With thermal switch and power fuse		

#### Material in contact with the pumped media

Type/OrderNo.	Pump head	Diaphragm	Valves
N 811 KN.18	PPS	EPDM	FPM
N 811 KT.18	PPS	PTFE-coated	FFPM

Motors with other voltages and frequencies on request.

<sup>1)</sup> at atm. pressure

## Accessories and spare parts


Description	Details	Order No.
Silencer	G 1/8	000345
Filter	G 1/8	000346
Hose connector	G 1/8, PVDF, for tube ID 6 mm	014052
Fine control valve with pressure gauge	pressure side	001786
Fine control valve with vacuum gauge	suction side	001787
Spares kit	for N 811 KN.18	044066
Spares kit	for N 811 KT.18	044067




# Diaphragm Pumps for Air, Gases and Vapours

## Dimensions and performance characteristics


Pump down time for 5 l receiver

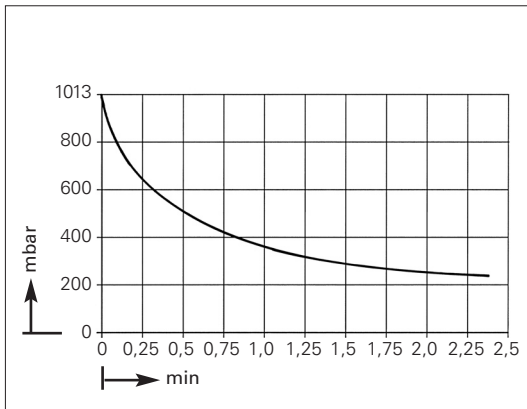

Performance characteristics

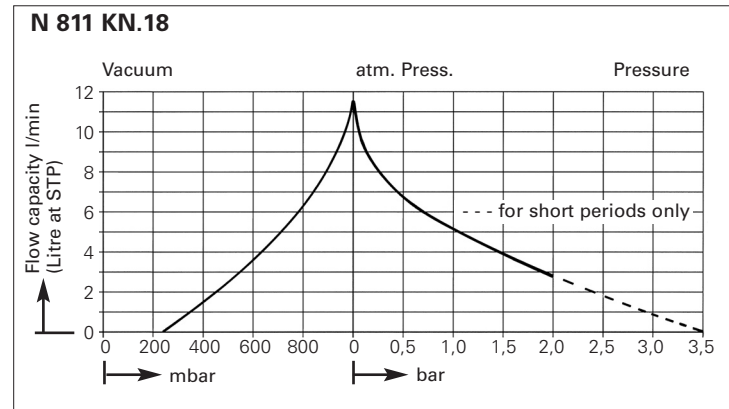

Performance characteristics

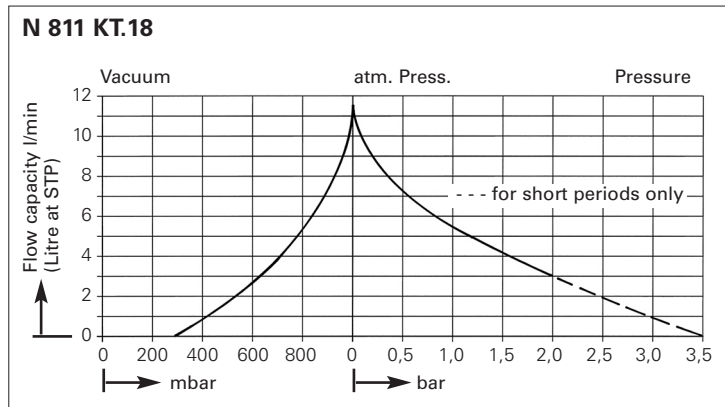

## Dimensions (mm)

